

LIBERTY

4-Channel, 12-Band Digital Amplifier 100% Digital Processing

Features and available options:

- Adaptive Feedback Rejection for Improved Stability, Higher Gains, and Larger Venting
- New Layered Noise Reduction Reduces Noise When Speech is not present
- 4-Channel WDRC Compression
- Excellent Performance in Noisy Environments
- Selectable Time Constant Modes
- Look-ahead Detection for Loud Impulse Noises
- Independent Expansion in Each Channel
- Optional Auto T-coil or Manual T-coil
- Can be either a regular BTE or Open fit BTE
- Available in 2 different color cases
Beige, and Beige with Cocoa Brown insert.
- 12-Band Equalizer for Fine Tuning at all Frequencies
- Up to four Programmable Memories with Memory Identification Tones
- Selection of Directional Patterns
- Wideband Output Compression limiting
- MPO Limiting Threshold Adjustments for each Channel
- Low Battery Warning Tone Alerts User of Low Battery Level
- A variety of tube and dome sizes available


	Ear Hook	Thin Tube
Battery Size	13A	13A
SSPL90 Max. (db SPL)	130	120
HF Average OSPL90 (db SPL)	125	110
Peak Gain (db)	60	35
HF Average Full on Gain (db)	55	30
Frequency Range (Hz)	200-8000	200-8000
Batter Life (Hours)	278	278

Electronic Digital Processing Option:

- Programmable Directional Patterns --
- Omni Directional Cardioid
- Super-cardioid Hyper-cardioid



**LIBERTY CONVERTIBLE
OTE/BTE**

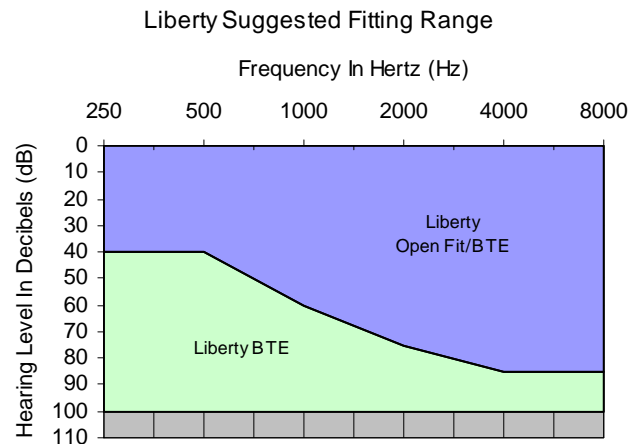

Highest Quality at Lowest Prices

Please Contact a Prairie Labs Sales Representative for more details (800) 322-8238